


# JumboMix<sup>®</sup> 3500

Large capacity lab blender

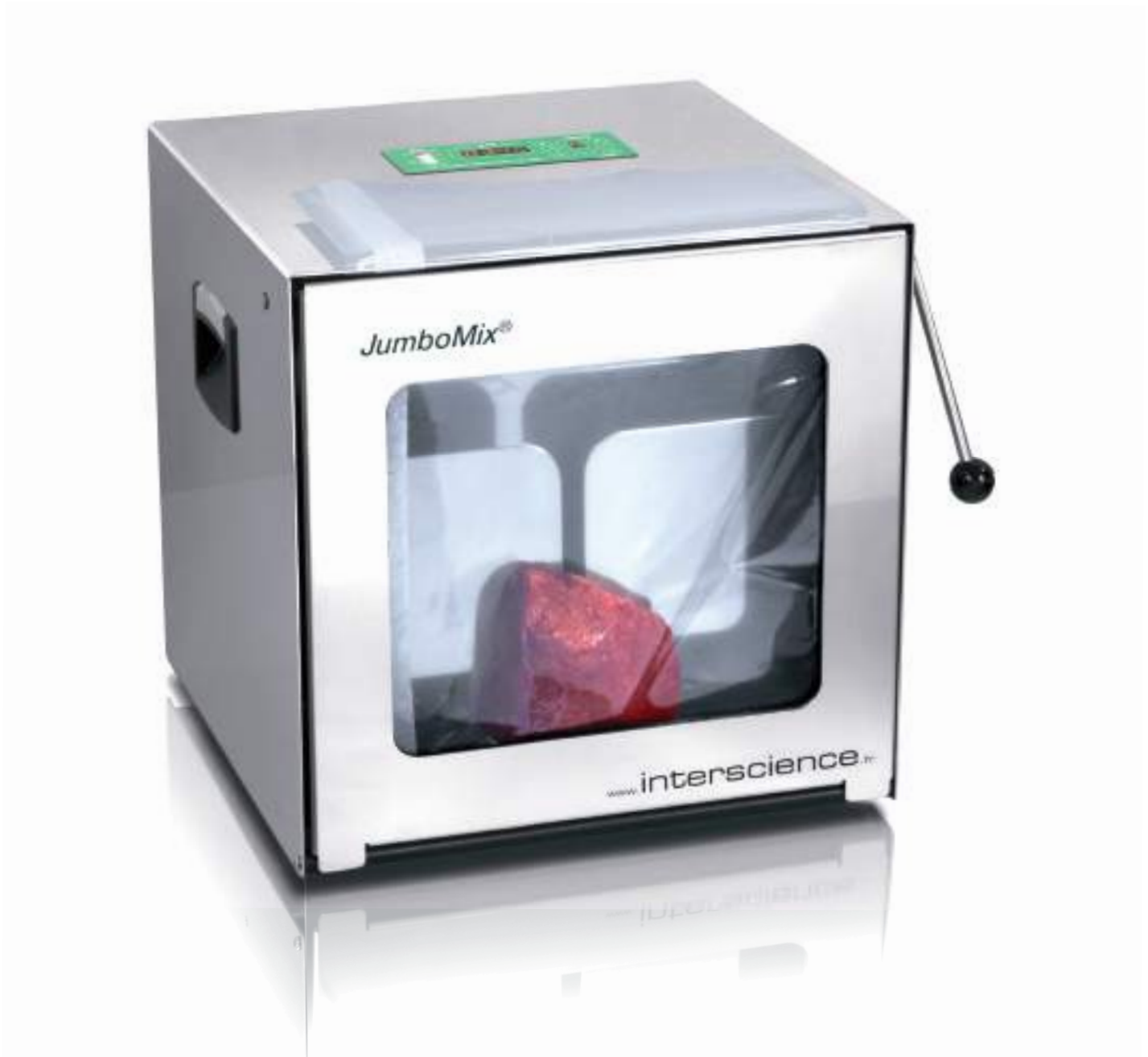

Accuracy and heavy duty blending  
for large samples

interscience

# Blend up to 5000 mL!

**JumboMix®** lab blender guarantees optimal bacterial extraction by the action of its 2 paddles. Adapted to any volume between 200 mL and 5000 mL, it guarantees the homogenization in 30 to 60 s of all samples.

**Applications:** E. Coli O157:H7, trichinella (WarmMix® model), salmonella, baby milk powder, chocolate, sterility control tests, Cryptosporidium and Giardia parasites, sample pooling (e.g.: 15 samples of 25 g to reach 375 g)...

## Powerful and heavy duty

- Up to 5000 mL capacity
- Ideal for your 375 g samples diluted to 1/10<sup>th</sup>
- Very silent

## Ergonomic

- Side-by-side paddle stop: easy insertion of the bag
- All 304L stainless steel trolley, ergonomic height (in option)
- Digital screen

## Easy cleaning

- All stainless steel
- 270° fully opening door

## Security

- The handle shuts and starts the blending
- **Q-Tight®** closing: no leaks upwards
- Security drip tray



## Guarantee

- Life-time guarantee: Shock absorbers and window
- 3 year guarantee (after validation of the card)
- Made in France

## PATENTED


### Adjustable blending power

- Adapt the blending power to your sample
- Up to 22.6 kg thrust: hard samples are blended




### Ideal for your applications

- 375 g for salmonella
- 75, 325 or 375 g sample for E. Coli O157:H7
- 125 g sample for Listeria and Listeria monocytogenes
- Sample pooling

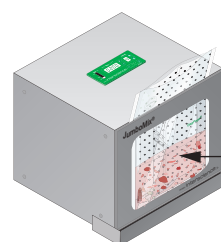

#### SAMPLE POOLING

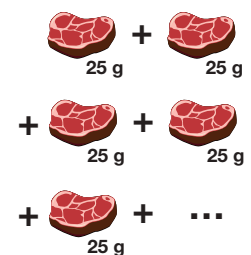

# JumboMix® 3500

3 models available:



## JumboMix® 3500 VP

> Essential features

Ref. 031 230



## JumboMix® 3500 VW

> Window door

Ref. 032 230



## JumboMix® 3500 WarmMix®

> Thermal regulation up to + 50°C

Ref. 033 230



Technical specifications at the back of the brochure →


# JumboLine®

Ideal for your preparation of 375 g samples diluted to 1/10<sup>th</sup>!



## JumboGravimat®

Gravimetric dilutor  
Ref. 512 235

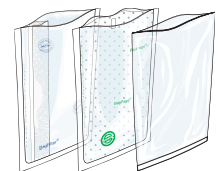

## BagFilter® 3500

Sterile bag with lateral filter  
Ref. 113 510

## BagPage® 3500

Sterile bag with full-page filter  
Ref. 123 010

## BagLight® 3500

Sterile bag  
Ref. 133 025



## JumboStand

Ergonomic trolley  
Ref. 037 000



## BagOpen® 3500

Bag opener  
Ref. 211 350



## BagRack® 3500

2 comp. - Ref. 221 352  
4 comp. - Ref. 221 350



## BagClip® 3500






Closing clips  
Ref. 231 350




## BagPipet® & BagTips

Pipet - Ref. 251 091  
Sterile tips - Ref. 252 024

# Technical specifications

	JumboMix® 3500 VP	JumboMix® 3500 VW	JumboMix® 3500 WarmMix®	
Reference	031 230	032 230	033 230	
<b>Equipment</b>	Window door	-	✓	-
	270° opening door	✓	✓	✓
	Security drip tray	in option	✓	✓
	Adjustable blending power (up to 22.6 kg thrust)	✓	✓	✓
	Blending power indicator	✓	✓	✓
	Digital screen	✓	✓	✓
	Q-Tight® anti-leak system	✓	✓	✓
	Side-by-side paddle stop	✓	✓	✓
	All stainless steel	✓	✓	✓
	Compact and ergonomic	✓	✓	✓
Regulated heating door	-	-	✓	
<b>Technical specifications</b>	Door weight	6.9 kg	3.9 kg	6.1 kg
	Blending time	from 1 s to 59 min, ∞		
	Speed	variable 1.5 / 3 / 4.5 / 6 strokes/s		
	Dimensions (w x d x h)	52 x 47 x 47 cm		63 x 50 x 47 cm
	Weight	58 kg	56 kg	63 kg
	Power	100-127 V~ / 220-240 V~ 50-60 Hz		
Capacity	200 - 5000 mL			
<b>Guarantee</b>	Made in France	✓	✓	✓
	Life guarantee	Shock absorbers and window door		
	3 years guarantee (after registration)	✓	✓	✓
<b>Certified production</b>	    			

## Discover our complete range for microbiology

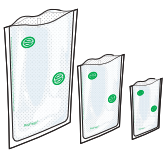

**BagFilter®**  
Sterile bags



**Gravimat®**  
Gravimetric dilutors



**FlexiPump®**  
Peristaltic pumps



**easySpiral Dilute®**  
Serial dilutors and platers



**Scan®**  
Colony counters

Distributor

### interscience INTERNATIONAL

30, ch. Bois Arpents - 78860 St Nom - FRANCE  
T: +33 (0)1 34 62 62 61 - F: +33 (0)1 34 62 43 03  
info@interscience.fr - www.interscience.fr

### interscience USA & CANADA

199 Weymouth ST. - ROCKLAND - MA 02370 - USA  
T: +1 781 792 2133 - F: +1 781 792 2134  
info@intersciencelab.com - www.intersciencelab.com

### interscience CHINA

上海市徐汇区吴兴路277号锦都大厦718室-200030  
电话: +86 (0)21-64739390 - 手机: +86 13040616422  
sales.china@interscience.fr - www.interscience.fr/china